

## Scheme 2024, Assessment Run 72

Aalborg, 2<sup>nd</sup> September 2024

Dear participant

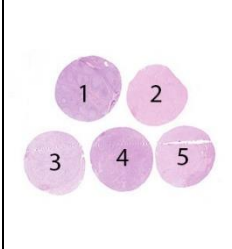

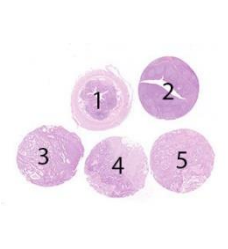
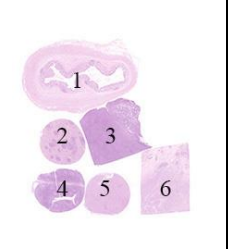
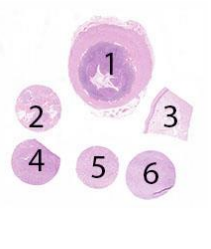
Pairs of unstained slides of NordiQC multi-tissue samples have been sent according to your submitted protocols for this run.

If you are in any doubt, please read the details on the website: [www.nordiqc.org](http://www.nordiqc.org), or send us an e-mail: [nordiqc@rn.dk](mailto:nordiqc@rn.dk). Also inform us if slides are missing.

Please note your participant number (NQC...) shown on the address label. This number **must** be indicated on all slides and used for any communication with NordiQC.

### Contents of the NordiQC multi-tissue samples

Top row closest to the glass label.

				
<b>Ki67:</b> 1: Tonsil 2: Liver 3: Seminoma 4: Neuroendocrine carcinoma, NOS 5: Pancreatic neuroendocrine tumor	<b>P16:</b> 1: Tonsil 2: Uterine Cervix 3: Cervical CIN3 4: Cervical squamous cell carcinoma 5-6: Head and neck squamous cell carcinoma	<b>MSH6:</b> 1: Appendix 2: Tonsil 3: Endometrial adenocarcinoma 4-5: Colon adenocarcinoma	<b>CD8:</b> 1: Appendix 2: Spleen 3: B-CLL 4: Tonsil 5: T-cell lymphoma 6: Breast carcinoma	<b>TRPS1:</b> 1: Appendix 2: Breast hyperplasia 3: Uterine cervix 4: NSCLC 5-6: TNBC

For each marker stain the sections according to the protocols you have uploaded on our website.

**If any changes are made in the protocols, upload the used protocols as these are used for data analysis and individual feed-back in case needed.**

In-house controls shall, if possible, be placed on the NordiQC slide (see next page).

Mark all slides with your NordiQC participant number.

For each marker return only one stained slide. Keep the remaining NordiQC slides to serve as own control set.

Do **not** return in-house controls on separate slides.

Deadline for slide return: **10<sup>th</sup> October**

Slides return address:  
 NordiQC  
 Institute of Pathology  
 Aalborg University Hospital  
 Ladegaardsgade 3  
 9000 Aalborg  
 Denmark

Overall and individual results are available on the website on 10<sup>th</sup> December.

Thank you for your assistance.

Best regards  
 NordiQC

General information;

## NordiQC tissue material and sections:

All included tissues in the NordiQC multi-tissue samples have been fixed in 10% neutral buffered formalin and embedded in paraffin.

The sections of the multi-tissue blocks are sectioned at 3-4  $\mu\text{m}$ .

After sectioning, the slides are air dried at room temp. for 48 hours, after which they are stored at  $-80^{\circ}\text{C}$  until distributed.

After receipt of slides, store them at  $-20^{\circ}\text{C}$  until staining.

The slides are unbaked.

Before IHC staining in-house control material should be mounted on the NordiQC slide and dried at e.g.  $37^{\circ}\text{C}$  overnight or  $55-60^{\circ}\text{C}$  for one hour to secure section adhesion.

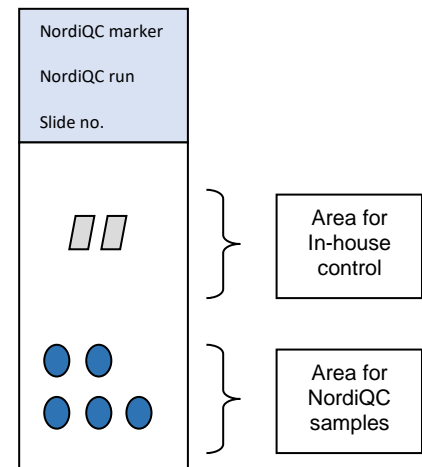


Fig. 1 Example of NordiQC slide

## Slide return

Return slides to NordiQC after IHC staining has been completed and in due time for deadlines indicated in this accompany letter and at [www.nordiqc.org/modules](http://www.nordiqc.org/modules). The slides must be returned in the same or a similar box in which it was received.

You can see if NordiQC has received your slides. Please log in at our website and choose "Protocol submission". When the slides are received at NordiQC and registered, this will be marked with the date received in the column "Slide received by NordiQC".

## Submitted sections

Submitted sections are stored by NordiQC for future documentation and will not be returned to the laboratories. However, laboratories can request their stained sections for review e.g. in case their own control set is lost.

## In-house on-slide controls

In-house on-slide controls shall serve as participants internal IHC system control to verify successful test. In-house on-slide controls will not be evaluated by NordiQC as expected patterns and content of tissue are not known and validated by NordiQC. The results in the controls will not affect the final individual assessment mark.

## Results

The overall assessment results for all markers are published in reports available at [www.nordiqc.org/epitope](http://www.nordiqc.org/epitope) at the time indicated at <https://nordiqc.org/modules.php>.

The assessment reports will provide an overview of assessment criteria, protocols, reagents and instruments used by the participants and results obtained. The report will provide guidelines for recommendable control materials and descriptions / illustrations of the expected staining patterns in these.

Protocols giving optimal results will be published on the NordiQC website. Best practice protocols will typically and if possible be listed for the main IHC systems used by the participants.

In case of an insufficient result assessed as borderline or poor, the participant will receive individual suggestions for protocol improvement.